

## Cy3 Goat Anti-Mouse IgG (H+L) Antibody

### 概述

产品名称	Cy3 Goat Anti-Mouse IgG (H+L) Antibody
宿主	Goat
靶标种属	Mouse
经测试应用	适用于: ICC/IF, IHC-Fr, IHC-P, ICC, Flow Cyt, ELISA
偶联物	Cy3. Ex: 552 nm, Em: 565 nm

### 性能


产品形式	Liquid
存储溶液	Preservative: 0.02% Sodium azide Constituents: 23% Glycerol, PBS, 1% BSA
存储说明	Shipped at 4°C. Store at 4°C short term (2 weeks). Upon delivery aliquot. Store at -20°C. Avoid freeze / thaw cycle. Stable for 12 months at -20°C. <b>Keep away from Light.</b>
浓度	1 mg/mL
纯度	Immunogen affinity purified
纯化说明	This antibody was isolated by affinity chromatography using antigen coupled to agarose beads.
克隆	Polyclonal
同种型	IgG

### 应用

应用	说明
ICC/IF	1/500-1/2000
IHC-Fr	Use at an assay dependent dilution.
IHC-P	1/100-1/500
Flow Cyt	1/250-1/1000
ELISA	Use at an assay dependent dilution.

### 注意事项

本品适合做单标实验。要多标 (multiple-labeling) 实验, 建议客户自己测试, 或使用与相近物种血清蛋白或者免疫球蛋白预先亲和吸附处理过的二抗。

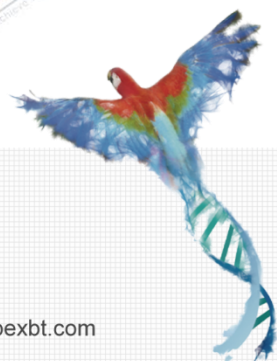

**APEx BIO Technology**

**[www.apexbt.com](http://www.apexbt.com)**

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: [info@apexbt.com](mailto:info@apexbt.com)