

Cy5 Goat Anti-Mouse IgG (H+L) Antibody

概述

产品名称	Cy5 Goat Anti-Mouse IgG (H+L) Antibody
宿主	Goat
靶标种属	Mouse
经测试应用	适用于: ICC/IF, IHC-Fr, IHC-P, ICC, Flow Cyt, ELISA
偶联物	Cy5. Ex: 650 nm, Em: 667 nm

性能


产品形式	Liquid
存储溶液	Preservative: 0.02% Sodium azide Constituents: 23% Glycerol, PBS, 1% BSA
存储说明	Shipped at 4°C. Store at 4°C short term (2 weeks). Upon delivery aliquot. Store at -20°C. Avoid freeze / thaw cycle. Stable for 12 months at -20°C. Keep away from Light.
浓度	1 mg/mL
纯度	Immunogen affinity purified
纯化说明	This antibody was isolated by affinity chromatography using antigen coupled to agarose beads.
克隆	Polyclonal
同种型	IgG

应用

应用	说明
ICC/IF	1/500-1/2000
IHC-Fr	Use at an assay dependent dilution.
IHC-P	1/100-1/500
Flow Cyt	1/250-1/1000
ELISA	Use at an assay dependent dilution.

注意事项

本品适合做单标实验。要多标 (multiple-labeling) 实验, 建议客户自己测试, 或使用与相近物种血清蛋白或者免疫球蛋白预先亲和吸附处理过的二抗。

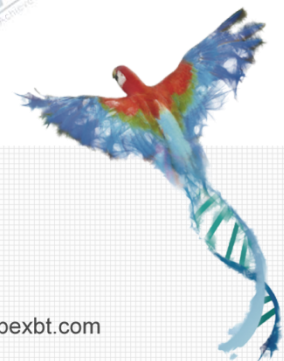

APExBIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com