

High Performance Liquid Chromatography

Analysis Information

Product Name: Recombinant Human ACE2 Protein
(Fc tag), HPLC-verified

Operator: David

Injection Date: 11/29/2022 12:16:14 AM

Batch No.: 1

HPLC Condition

Column : Athena C18, 3 μ m, 2.1mm \times 100 mm

Solvent A: ACN+0.1%TFA

Solvent B: H2O+0.1%TFA

Gradient:


Time	A	B
0.0min	0%	100%
1.5min	0%	100%
6.0min	95%	5%
9.0min	95%	5%
9.0min	Stop	Stop

Flow rate: 0.46ml/min

Wavelength: 214nm

Volume: 3 μ l

Result

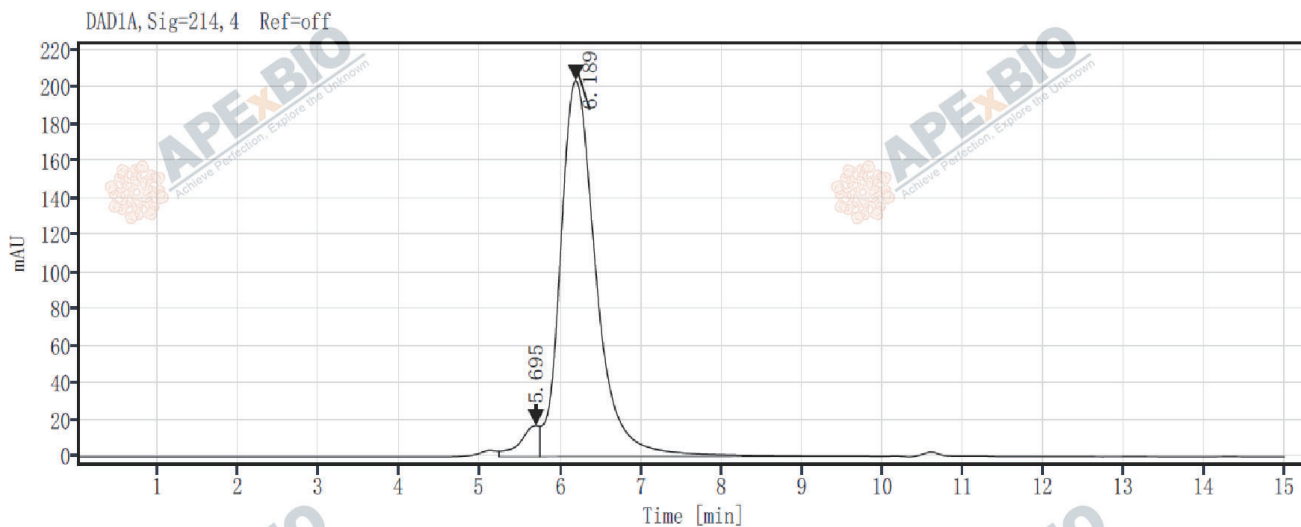

Rank	RetTime(min)	Area(mAU*s)	Area %
1	5.695	275.95	4.21
2	6.189	6283.76	95.79

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APEX BIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com